



Patient name Phantom Lunge
Patient ID 012345678910
Treatment plan name 10x10_9

Plan last save time 13 Aug 2015, 09:21:49 (hr:min:sec)
Report creation time 13 Aug 2015, 16:35:24 (hr:min:sec)
Plan and structure set approved No
Plan approved by -
Plan approval time -

Plan Report

Patient data

Patient ID	012345678910
Patient name	Phantom Lunge
Patient gender	Other
Patient birth date	08 Dec 1987
Treatment planning system	RayStation 4.5.0.19
Structure set UID	1.2.826.0.1.3680043.8.176.201581392252434.246.5555341432
Structure set approval data	
Approved	No
Approved by	-
Approval time	-

Treatment plan data

Treatment plan name	10x10_9
Plan last save time	13 Aug 2015, 09:21:49 (hr:min:sec)
Planned by	
Number of beam sets	1
Patient treatment position	HFS : Head First Supine
Treatment plan approval data	
Approved	No
Approved by	-
Approval time	-
Plan comment	
Planning image set	CT 62
CT to density table	DKFZ_HLUT 03 Apr 2012, 10:28:06 (hr:min:sec)
Patient scanning position	HFS
External ROI	External

General data

Treatment planning system	RayStation 4.5.0.19
Report creation time	13 Aug 2015, 16:35:24 (hr:min:sec)
Template name	RayStation treatment plan report
Patient coordinate system	IEC 61217

Density override

No density override

Beam Set overview

Beam Set name	10x10_9
Treatment technique	3D-CRT
Treatment unit	ARTISTE3
Number of beams	9

Warnings [10x10_9]

- The geometry 'ITV2_3_4' is derived and depends on geometries that are empty.
The geometry 'PTV_CT38' is derived and depends on geometries that are empty.
The geometry 'PTV_38' is derived and depends on geometries that are empty.
The geometry 'ITV_Ct50' is derived and depends on geometries that are empty.
The geometry 'PTV_50' is derived and depends on geometries that are empty.

Signatures

Signature 1 (Name/Signature/Date)

Signature 2 (Name/Signature/Date)

Beam Set Report

Beam Set data

Beam Set name	10x10_9
Modality	Photons
Treatment technique	3D-CRT
Number of beams	9
Number of segments	9
DICOM Plan UID	1.2.826.0.1.3680043.8.176.201581392149280.234.1872358020
Planning image set	CT 62
CT to density table	DKFZ_HLUT 03 Apr 2012, 10:28:06 (hr:min:sec)
Treatment unit	ARTISTE3
Commission time	05 Nov 2014, 10:34:16 (hr:min:sec)
Treatment machine scale	IEC 61217
Jaw labeling standard	IEC 61217
Energy [MV]	6.00
Dose calculation algorithm	Collapsed Cone, Version 3.0
Density calculation algorithm version	2.0
MU per fraction	525.25
Number of fractions	1
ROI(s) with density override	
Beam set approval data	
Approved	No
Approved by	-
Approval time	-

Beam Data Overview [Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99]

#	Beam name (Description)	Number of segments	Maximum jaw aperture [cm]		Gantry angle [deg]	Coll. angle [deg]	Couch angle [deg]	MU per fraction	Bolus [Y/N]	Block [Y/N]
			Y1	Y2						
1	1-B0 (1-B0)	1	-5.00	5.00	0.0	0.0	0.0	58.36	N	N
2	1-B1 (1-B1)	1	-5.00	5.00	40.0	0.0	0.0	58.36	N	N
3	1-B2 (1-B2)	1	-5.00	5.00	80.0	0.0	0.0	58.36	N	N
4	1-B3 (1-B3)	1	-5.00	5.00	120.0	0.0	0.0	58.36	N	N
5	1-B4 (1-B4)	1	-5.00	5.00	160.0	0.0	0.0	58.36	N	N
6	1-B5 (1-B5)	1	-5.00	5.00	200.0	0.0	0.0	58.36	N	N
7	1-B6 (1-B6)	1	-5.00	5.00	240.0	0.0	0.0	58.36	N	N
8	1-B7 (1-B7)	1	-5.00	5.00	280.0	0.0	0.0	58.36	N	N
9	1-B8 (1-B8)	1	-5.00	5.00	320.0	0.0	0.0	58.36	N	N

Objectives

No objectives defined

Constraints

No constraints defined

Prescription

Prescription	4.00 Gy to point dose at ●iso
Value [Gy]	4.00
Fulfilled	●Yes
Relates to beam set dose	

Patient setup



Localization point	
Treatment position	HFS : Head First Supine
POI	●RefPunkt
Position [cm]	X(Right-Left) = 0.15 , Y(Inf-Sup) = 3.07 , Z(Post-Ant) = 17.1
Patient setup	
Beams	1-B0, 1-B1, 1-B2, 1-B3, 1-B4, 1-B5, 1-B6, 1-B7, 1-B8


Isocenter [cm] X(R-L) = 0.44 , Y(I-S) = 0.47 , Z(P-A) = 25.99
Localization point - Isocenter [cm] X(R-L) = -0.29 , Y(I-S) = 2.6 , Z(P-A) = -8.89

Position patient such that lasers line up with patient marks.
Move the couch according to the PATIENT coordinate system:
RIGHT 0.29 cm (patient's right)
SUPERIOR 2.6 cm
POSTERIOR 8.89 cm

Beamset dose data

Isocenter [cm] Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99
Dose grid resolution [cm] Right-Left: 0.20 Inf-Sup: 0.20 Post-Ant: 0.20
Beams 1-B0, 1-B1, 1-B2, 1-B3, 1-B4, 1-B5, 1-B6, 1-B7, 1-B8

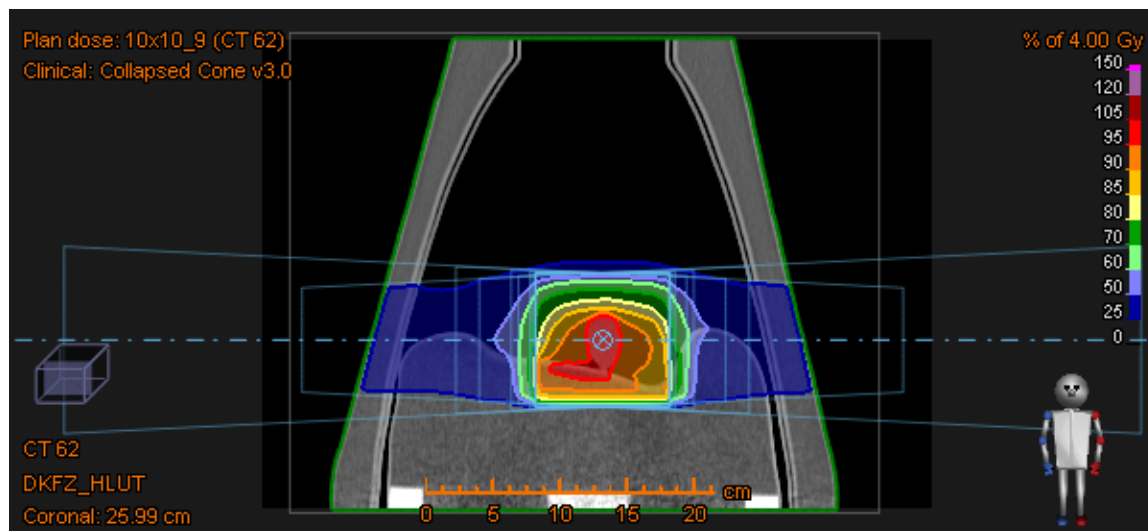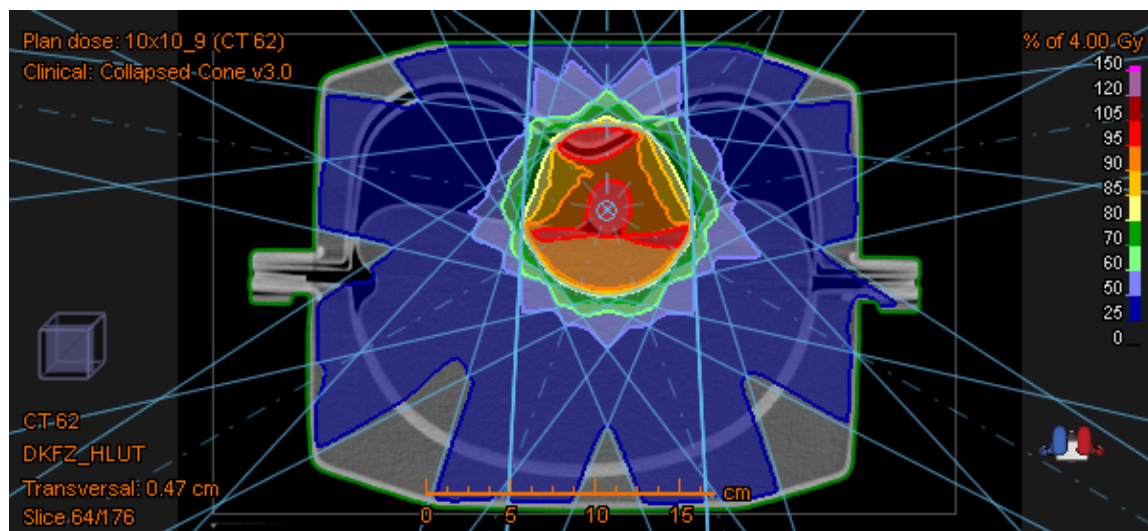

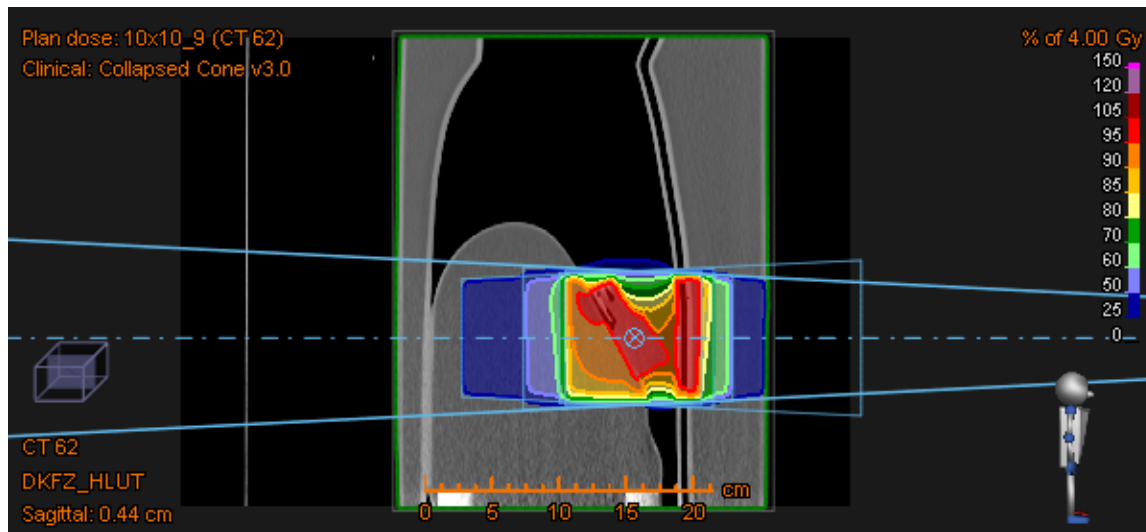


Points Of Interest

	Name	RefPunkt	Beam isocenters [cm]	Point - Isocenter [cm]
●	Type Dose [Gy] Location [cm]	Localization point 1.70 [Interpolated] Right-Left: 0.15 Inf-Sup: 3.07 Post-Ant: 17.1	Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99	Right-Left: -0.29 Inf-Sup: 2.60 Post-Ant: -8.89
●	Name Type Dose [Gy] Location [cm]	Iso Isocenter 0.00 [Interpolated] N/A		
●	Name Type Dose [Gy] Location [cm]	iso Isocenter 4.00 [Interpolated] Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99	Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99	Right-Left: 0.00 Inf-Sup: 0.00 Post-Ant: 0.00

Clinical goals

There are no clinical goals

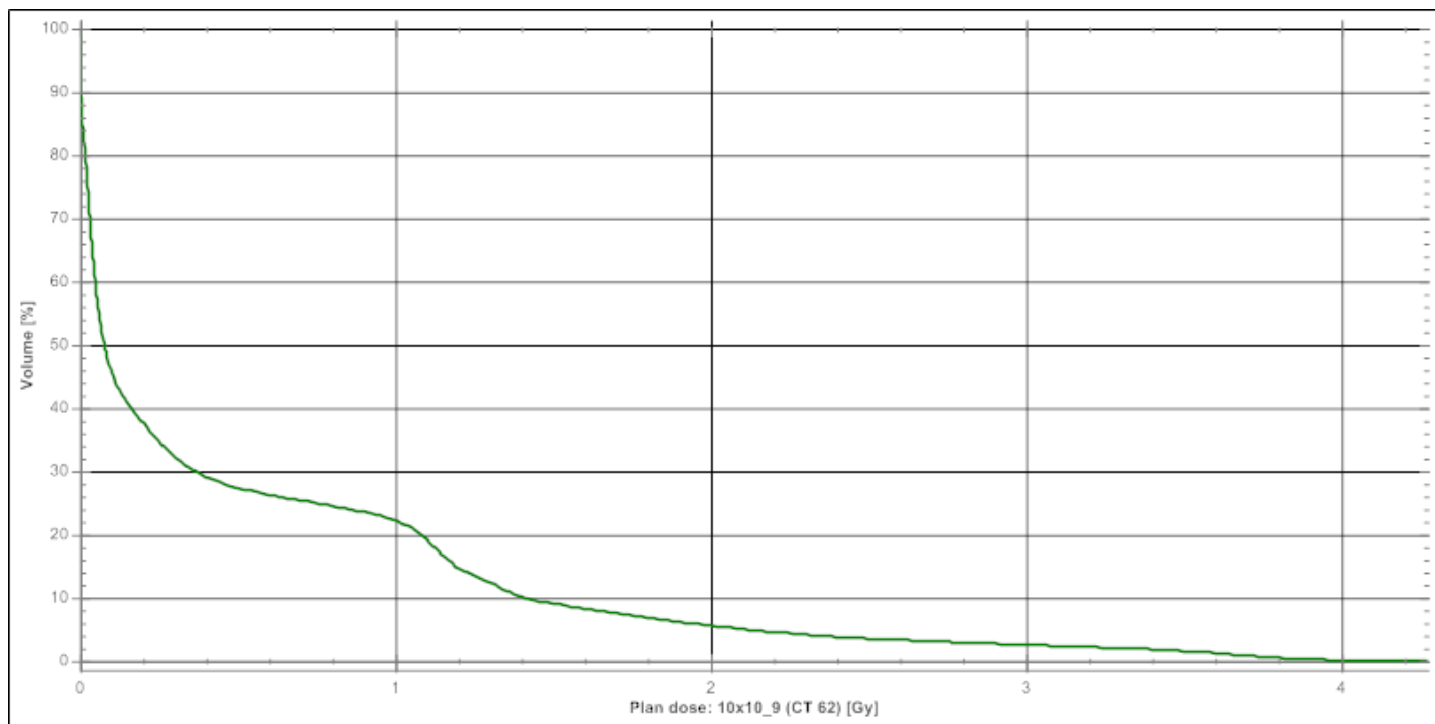

POI Dose statistics [Beam Set dose]

Dose	POI	Dose [Gy]	Position Right-Left: [cm]	Inf-Sup: [cm]	Post-Ant: [cm]
Plan dose: 10x10_9 (CT 62)	● RefPunkt	1.70	0.15	3.07	17.1
Plan dose: 10x10_9 (CT 62)	● Iso	-	-	-	-
Plan dose: 10x10_9 (CT 62)	● iso	4.00	0.44	0.47	25.99

ROI Dose statistics [Beam Set dose]

Name	Volume [cm³]	D99 [Gy]	D98 [Gy]	D95 [Gy]	Average [Gy]	D50 [Gy]	D2 [Gy]	D1 [Gy]	% outside grid
External	28676.68	0.00	0.00	0.00	0.48	0.08	3.32	3.66	0
GTV									-
ITV									-
ITV_38									-
ITV_50									-
ITV_CT38									-
ITV_Ct50									-
ITV2									-
ITV2_3									-
ITV2_3_4									-
ITV3									-
ITV4									-
PTV									-
PTV_38									-
PTV_50									-
PTV_CT 38->CT 41									-
PTV_CT 38->CT 42									-
PTV_CT 38->CT 43									-
PTV_CT 38->CT 44									-



Patient name

Patient ID

Treatment plan name

Phantom Lunge

012345678910

10x10_9

Plan last save time

Report creation time

Plan and structure set approved

Plan approved by

Plan approval time

13 Aug 2015, 09:21:49 (hr:min:sec)

13 Aug 2015, 16:35:24 (hr:min:sec)

No

-

-

PTV_CT 38->CT 44 (1)									-
PTV_CT 38->CT 45									-
PTV_CT 38->CT 46									-
PTV_CT 38->CT 46 (1)									-
PTV_CT 38->CT 47									-
PTV_CT 38->CT 48									-
PTV_CT 38->CT 49									-
PTV_CT 38->CT 49 (1)									-
PTV_CT38									-
Tumor									-

External



This ROI is set as the external ROI that defines the outer border of the patient

Beam data

Beam name	1-B0
Beam number	1
Beam description	1-B0
Patient coordinate system	IEC 61217
Isocenter [cm]	Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99
Gantry angle [deg]	0.0
Collimator angle [deg]	0.0
Couch angle [deg]	0.0
Treatment technique	3D-CRT
Number of fractions	1
Beam MU/fraction	58.36
Total beam MU	58.36
Beam weight	0.11
Number of segments	1
Dose calculation algorithm	Collapsed Cone, Version 3.0
Treatment unit	ARTISTE3
Commission time	05 Nov 2014, 10:34:16 (hr:min:sec)
Energy [MV]	6.00
Jaw max aperture width [cm]	-
X1 [cm]	-
X2 [cm]	-
Jaw max aperture height [cm]	10.00
Y1 [cm]	-5.00
Y2 [cm]	5.00
Source to skin distance (isocenter) [cm]	90.09
Source to surface distance (isocenter) [cm]	90.09
Bolus data	
No bolus	

Beam dose specification point

Coordinates [cm]	Isocenter
Dose per fraction [Gy]	0.501
Physical depth [cm]	9.91
Water equivalent depth [cm]	8.58
Source to skin distance [cm]	90.09
Source to surface distance [cm]	90.09

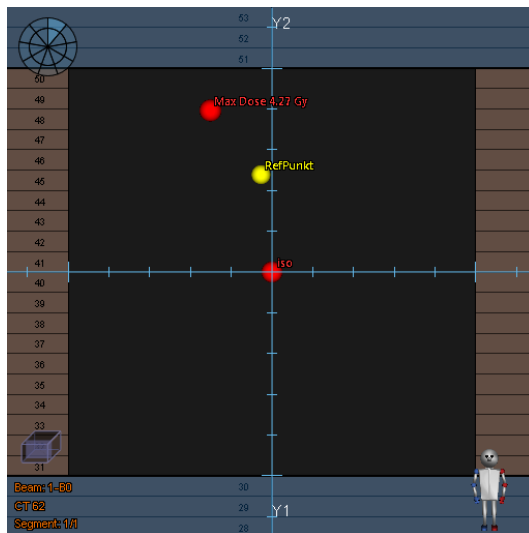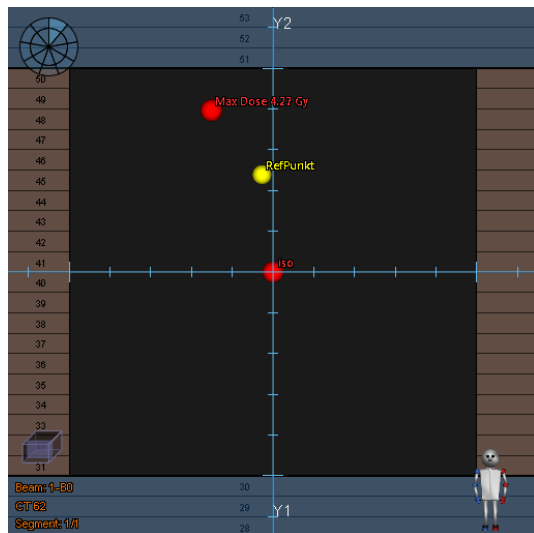

Segments



Seg. No.	MU/Fraction	Jaw positions [cm]	
		Y1	Y2
1	58.36	-5.00	5.00

Beam data

Beam name	1-B1
Beam number	2
Beam description	1-B1
Patient coordinate system	IEC 61217
Isocenter [cm]	Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99
Gantry angle [deg]	40.0
Collimator angle [deg]	0.0
Couch angle [deg]	0.0
Treatment technique	3D-CRT
Number of fractions	1
Beam MU/fraction	58.36
Total beam MU	58.36
Beam weight	0.11
Number of segments	1
Dose calculation algorithm	Collapsed Cone, Version 3.0
Treatment unit	ARTISTE3
Commission time	05 Nov 2014, 10:34:16 (hr:min:sec)
Energy [MV]	6.00
Jaw max aperture width [cm]	-
X1 [cm]	-
X2 [cm]	-
Jaw max aperture height [cm]	10.00
Y1 [cm]	-5.00
Y2 [cm]	5.00
Source to skin distance (isocenter) [cm]	87.36
Source to surface distance (isocenter) [cm]	87.36
Bolus data	
No bolus	

Beam dose specification point

Coordinates [cm]	Isocenter
Dose per fraction [Gy]	0.522
Physical depth [cm]	12.64
Water equivalent depth [cm]	6.86
Source to skin distance [cm]	87.36
Source to surface distance [cm]	87.36

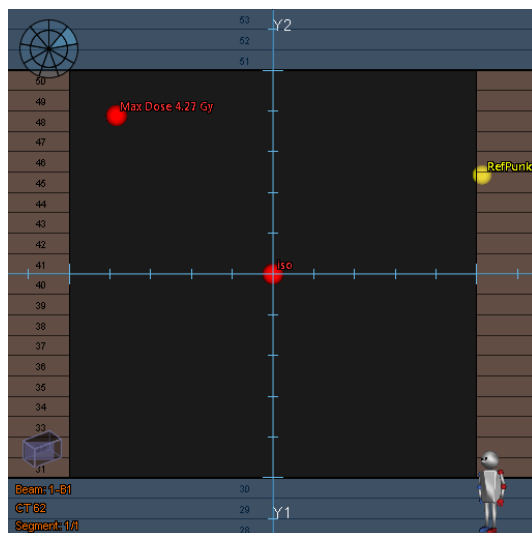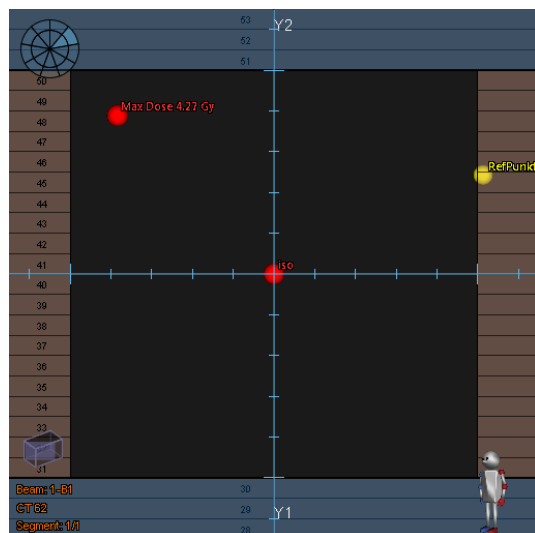

Segments



Seg. No.	MU/Fraction	Jaw positions [cm]	
		Y1	Y2
1	58.36	-5.00	5.00

Beam data

Beam name	1-B2
Beam number	3
Beam description	1-B2
Patient coordinate system	IEC 61217
Isocenter [cm]	Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99
Gantry angle [deg]	80.0
Collimator angle [deg]	0.0
Couch angle [deg]	0.0
Treatment technique	3D-CRT
Number of fractions	1
Beam MU/fraction	58.36
Total beam MU	58.36
Beam weight	0.11
Number of segments	1
Dose calculation algorithm	Collapsed Cone, Version 3.0
Treatment unit	ARTISTE3
Commission time	05 Nov 2014, 10:34:16 (hr:min:sec)
Energy [MV]	6.00
Jaw max aperture width [cm]	-
X1 [cm]	-
X2 [cm]	-
Jaw max aperture height [cm]	10.00
Y1 [cm]	-5.00
Y2 [cm]	5.00
Source to skin distance (isocenter) [cm]	84.72
Source to surface distance (isocenter) [cm]	84.72
Bolus data	
No bolus	

Beam dose specification point

Coordinates [cm]	Isocenter
Dose per fraction [Gy]	0.550
Physical depth [cm]	15.28
Water equivalent depth [cm]	5.42
Source to skin distance [cm]	84.72
Source to surface distance [cm]	84.72

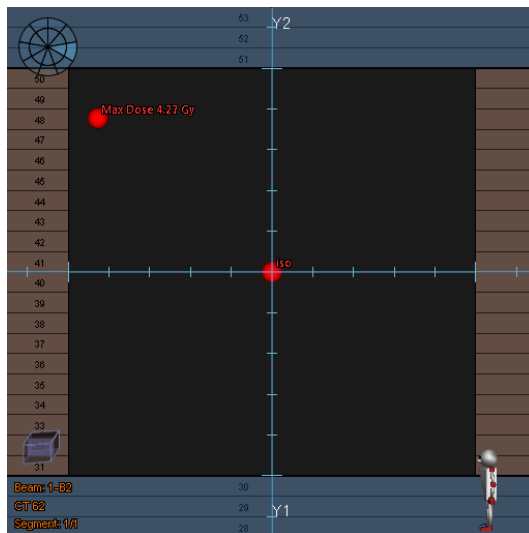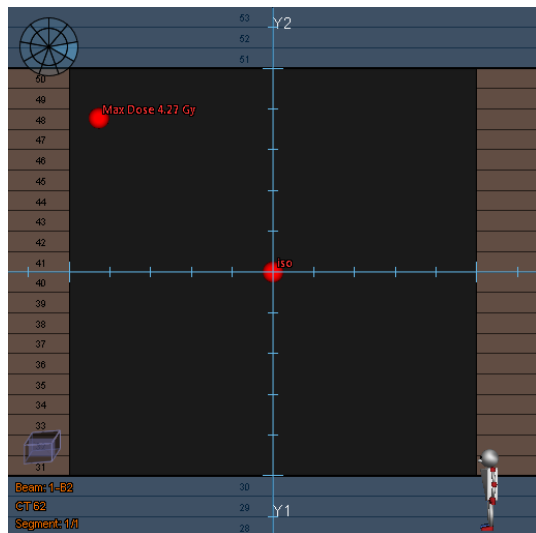

Segments



Seg. No.	MU/Fraction	Jaw positions [cm]	
		Y1	Y2
1	58.36	-5.00	5.00

Beam data

Beam name	1-B3
Beam number	4
Beam description	1-B3
Patient coordinate system	IEC 61217
Isocenter [cm]	Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99
Gantry angle [deg]	120.0
Collimator angle [deg]	0.0
Couch angle [deg]	0.0
Treatment technique	3D-CRT
Number of fractions	1
Beam MU/fraction	58.36
Total beam MU	58.36
Beam weight	0.11
Number of segments	1
Dose calculation algorithm	Collapsed Cone, Version 3.0
Treatment unit	ARTISTE3
Commission time	05 Nov 2014, 10:34:16 (hr:min:sec)
Energy [MV]	6.00
Jaw max aperture width [cm]	-
X1 [cm]	-
X2 [cm]	-
Jaw max aperture height [cm]	10.00
Y1 [cm]	-5.00
Y2 [cm]	5.00
Source to skin distance (isocenter) [cm]	82.81
Source to surface distance (isocenter) [cm]	82.81
Bolus data	
No bolus	

Beam dose specification point

Coordinates [cm]	Isocenter
Dose per fraction [Gy]	0.374
Physical depth [cm]	17.19
Water equivalent depth [cm]	16.04
Source to skin distance [cm]	82.81
Source to surface distance [cm]	82.81

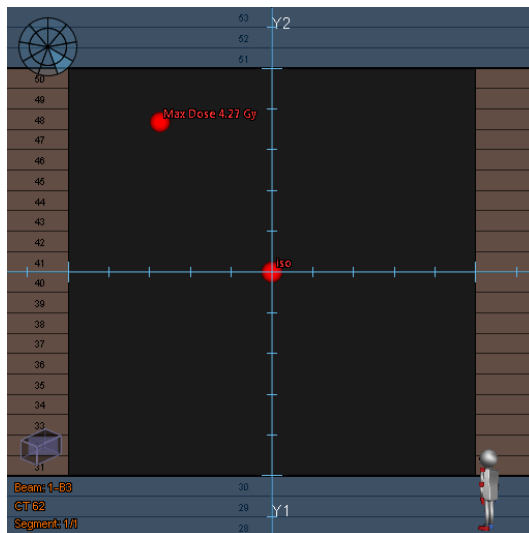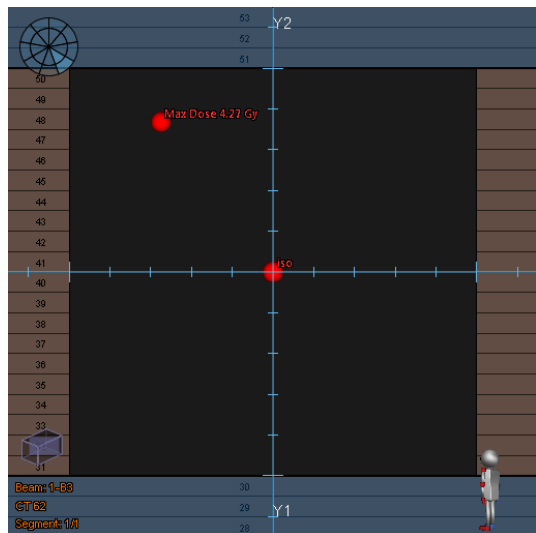

Segments



Seg. No.	MU/Fraction	Jaw positions [cm]	
		Y1	Y2
1	58.36	-5.00	5.00

Beam data

Beam name	1-B4
Beam number	5
Beam description	1-B4
Patient coordinate system	IEC 61217
Isocenter [cm]	Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99
Gantry angle [deg]	160.0
Collimator angle [deg]	0.0
Couch angle [deg]	0.0
Treatment technique	3D-CRT
Number of fractions	1
Beam MU/fraction	58.36
Total beam MU	58.36
Beam weight	0.11
Number of segments	1
Dose calculation algorithm	Collapsed Cone, Version 3.0
Treatment unit	ARTISTE3
Commission time	05 Nov 2014, 10:34:16 (hr:min:sec)
Energy [MV]	6.00
Jaw max aperture width [cm]	-
X1 [cm]	-
X2 [cm]	-
Jaw max aperture height [cm]	10.00
Y1 [cm]	-5.00
Y2 [cm]	5.00
Source to skin distance (isocenter) [cm]	81.18
Source to surface distance (isocenter) [cm]	81.18
Bolus data	No bolus

Beam dose specification point

Coordinates [cm]	Isocenter
Dose per fraction [Gy]	0.344
Physical depth [cm]	18.82
Water equivalent depth [cm]	18.76
Source to skin distance [cm]	81.18
Source to surface distance [cm]	81.18

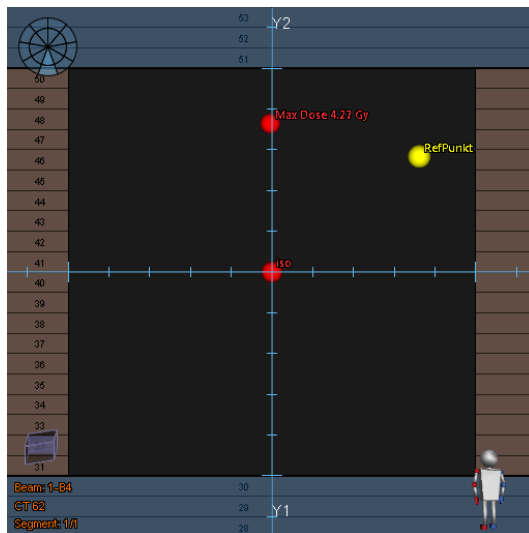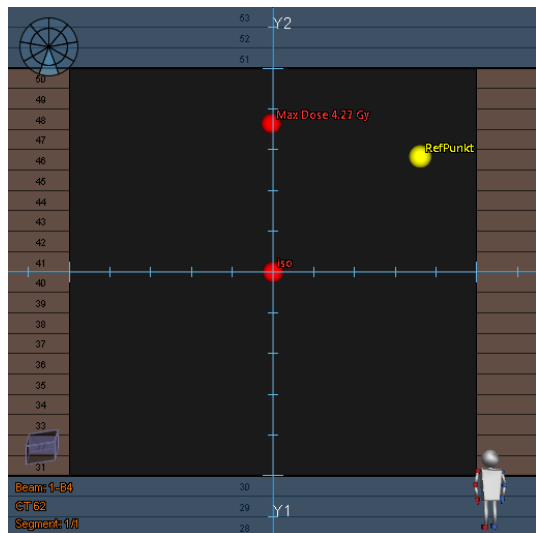

Segments



Seg. No.	MU/Fraction	Jaw positions [cm]	
		Y1	Y2
1	58.36	-5.00	5.00

Beam data

Beam name	1-B5
Beam number	6
Beam description	1-B5
Patient coordinate system	IEC 61217
Isocenter [cm]	Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99
Gantry angle [deg]	200.0
Collimator angle [deg]	0.0
Couch angle [deg]	0.0
Treatment technique	3D-CRT
Number of fractions	1
Beam MU/fraction	58.36
Total beam MU	58.36
Beam weight	0.11
Number of segments	1
Dose calculation algorithm	Collapsed Cone, Version 3.0
Treatment unit	ARTISTE3
Commission time	05 Nov 2014, 10:34:16 (hr:min:sec)
Energy [MV]	6.00
Jaw max aperture width [cm]	-
X1 [cm]	-
X2 [cm]	-
Jaw max aperture height [cm]	10.00
Y1 [cm]	-5.00
Y2 [cm]	5.00
Source to skin distance (isocenter) [cm]	81.19
Source to surface distance (isocenter) [cm]	81.19
Bolus data	
No bolus	

Beam dose specification point

Coordinates [cm]	Isocenter
Dose per fraction [Gy]	0.344
Physical depth [cm]	18.81
Water equivalent depth [cm]	18.83
Source to skin distance [cm]	81.19
Source to surface distance [cm]	81.19

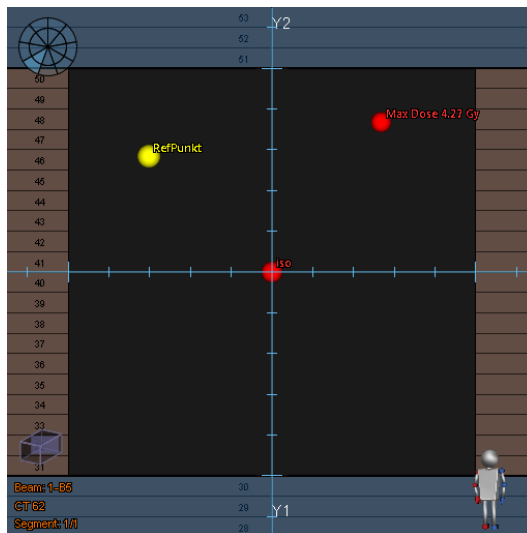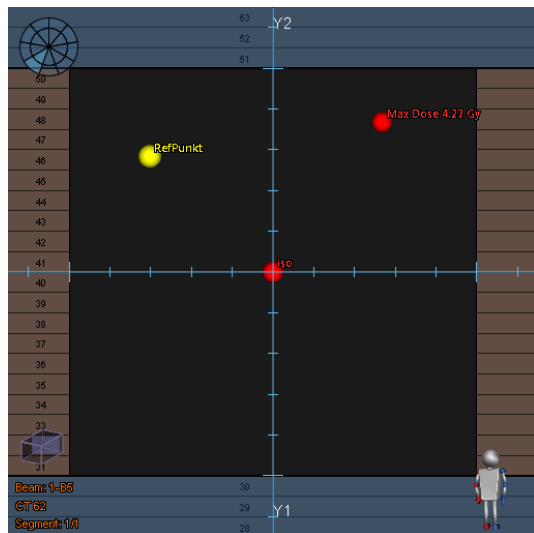

Segments



Seg. No.	MU/Fraction	Jaw positions [cm]	
		Y1	Y2
1	58.36	-5.00	5.00

Beam data

Beam name	1-B6
Beam number	7
Beam description	1-B6
Patient coordinate system	IEC 61217
Isocenter [cm]	Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99
Gantry angle [deg]	240.0
Collimator angle [deg]	0.0
Couch angle [deg]	0.0
Treatment technique	3D-CRT
Number of fractions	1
Beam MU/fraction	58.36
Total beam MU	58.36
Beam weight	0.11
Number of segments	1
Dose calculation algorithm	Collapsed Cone, Version 3.0
Treatment unit	ARTISTE3
Commission time	05 Nov 2014, 10:34:16 (hr:min:sec)
Energy [MV]	6.00
Jaw max aperture width [cm]	-
X1 [cm]	-
X2 [cm]	-
Jaw max aperture height [cm]	10.00
Y1 [cm]	-5.00
Y2 [cm]	5.00
Source to skin distance (isocenter) [cm]	79.51
Source to surface distance (isocenter) [cm]	79.51
Bolus data	
No bolus	

Beam dose specification point

Coordinates [cm]	Isocenter
Dose per fraction [Gy]	0.325
Physical depth [cm]	20.49
Water equivalent depth [cm]	19.56
Source to skin distance [cm]	79.51
Source to surface distance [cm]	79.51

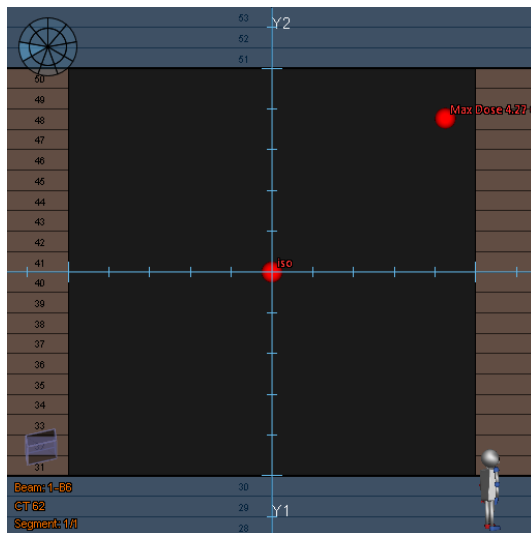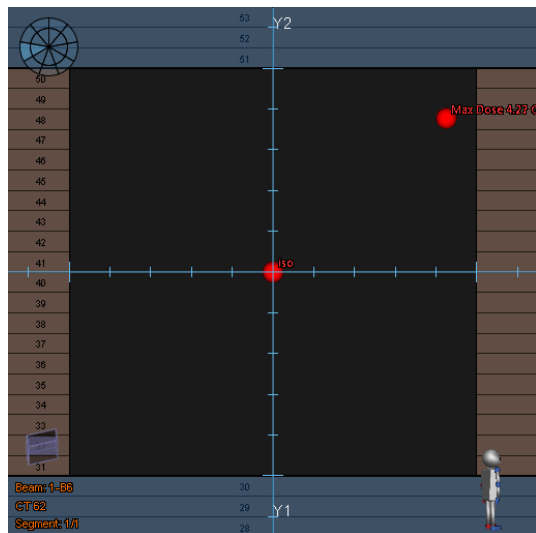

Segments


Seg. No.	MU/Fraction	Jaw positions [cm]	
		Y1	Y2
1	58.36	-5.00	5.00

Beam data

Beam name	1-B7
Beam number	8
Beam description	1-B7
Patient coordinate system	IEC 61217
Isocenter [cm]	Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99
Gantry angle [deg]	280.0
Collimator angle [deg]	0.0
Couch angle [deg]	0.0
Treatment technique	3D-CRT
Number of fractions	1
Beam MU/fraction	58.36
Total beam MU	58.36
Beam weight	0.11
Number of segments	1
Dose calculation algorithm	Collapsed Cone, Version 3.0
Treatment unit	ARTISTE3
Commission time	05 Nov 2014, 10:34:16 (hr:min:sec)
Energy [MV]	6.00
Jaw max aperture width [cm]	-
X1 [cm]	-
X2 [cm]	-
Jaw max aperture height [cm]	10.00
Y1 [cm]	-5.00
Y2 [cm]	5.00
Source to skin distance (isocenter) [cm]	82.45
Source to surface distance (isocenter) [cm]	82.45
Bolus data	
No bolus	

Beam dose specification point

Coordinates [cm]	Isocenter
Dose per fraction [Gy]	0.563
Physical depth [cm]	17.55
Water equivalent depth [cm]	4.61
Source to skin distance [cm]	82.45
Source to surface distance [cm]	82.45

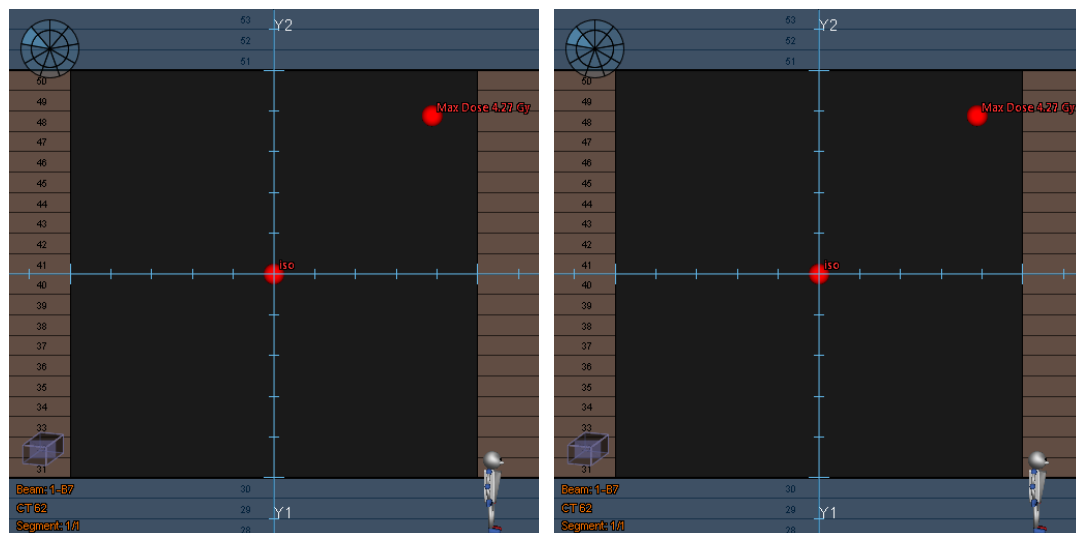

Segments



Seg. No.	MU/Fraction	Jaw positions [cm]	
		Y1	Y2
1	58.36	-5.00	5.00

Beam data

Beam name	1-B8
Beam number	9
Beam description	1-B8
Patient coordinate system	IEC 61217
Isocenter [cm]	Right-Left: 0.44 Inf-Sup: 0.47 Post-Ant: 25.99
Gantry angle [deg]	320.0
Collimator angle [deg]	0.0
Couch angle [deg]	0.0
Treatment technique	3D-CRT
Number of fractions	1
Beam MU/fraction	58.36
Total beam MU	58.36
Beam weight	0.11
Number of segments	1
Dose calculation algorithm	Collapsed Cone, Version 3.0
Treatment unit	ARTISTE3
Commission time	05 Nov 2014, 10:34:16 (hr:min:sec)
Energy [MV]	6.00
Jaw max aperture width [cm]	-
X1 [cm]	-
X2 [cm]	-
Jaw max aperture height [cm]	10.00
Y1 [cm]	-5.00
Y2 [cm]	5.00
Source to skin distance (isocenter) [cm]	87.09
Source to surface distance (isocenter) [cm]	87.09
Bolus data	
No bolus	




Beam dose specification point




Coordinates [cm]	Isocenter
Dose per fraction [Gy]	0.477
Physical depth [cm]	12.91
Water equivalent depth [cm]	9.06
Source to skin distance [cm]	87.09
Source to surface distance [cm]	87.09


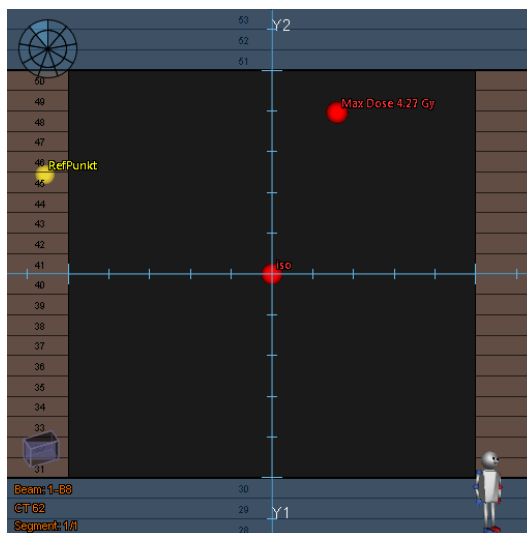
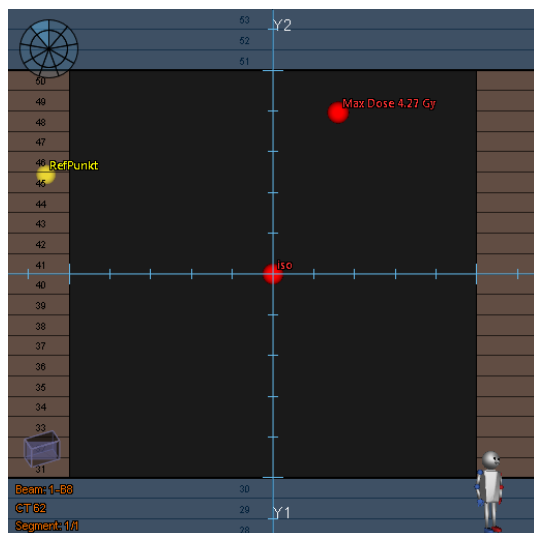


Segments

Seg. No.	MU/Fraction	Jaw positions [cm]	
		Y1	Y2
1	58.36	-5.00	5.00

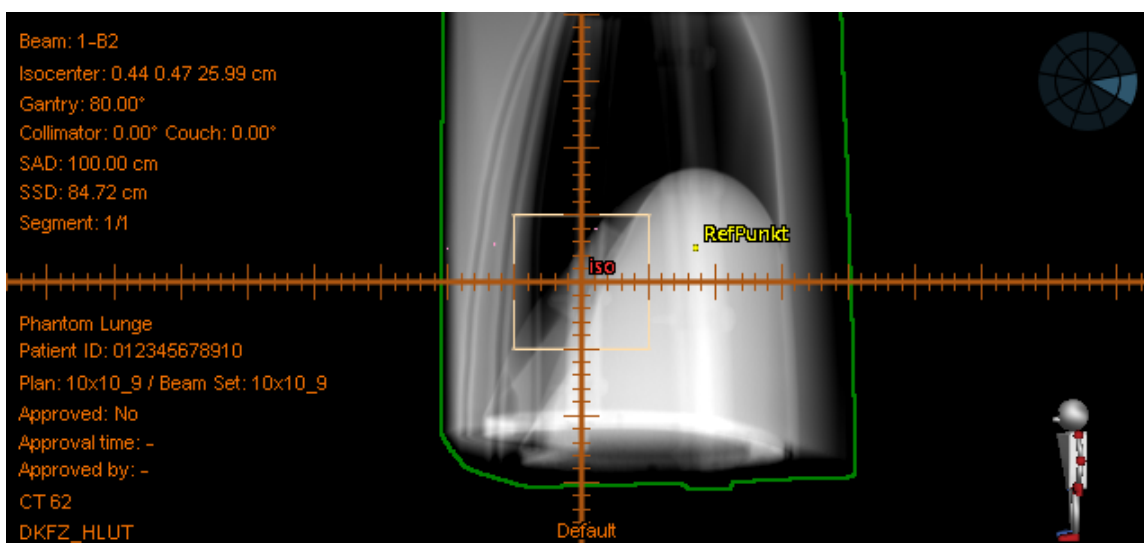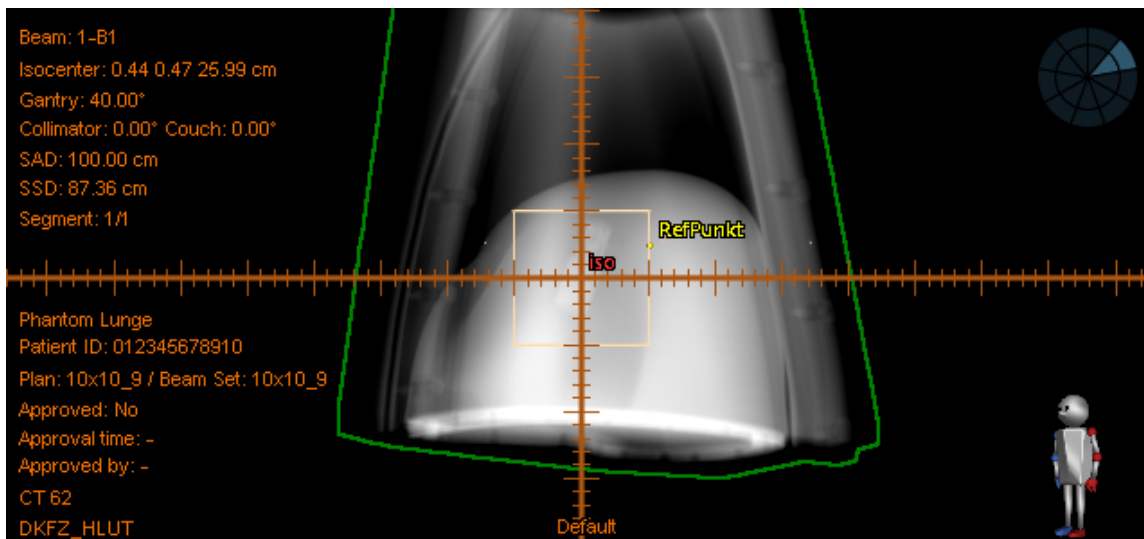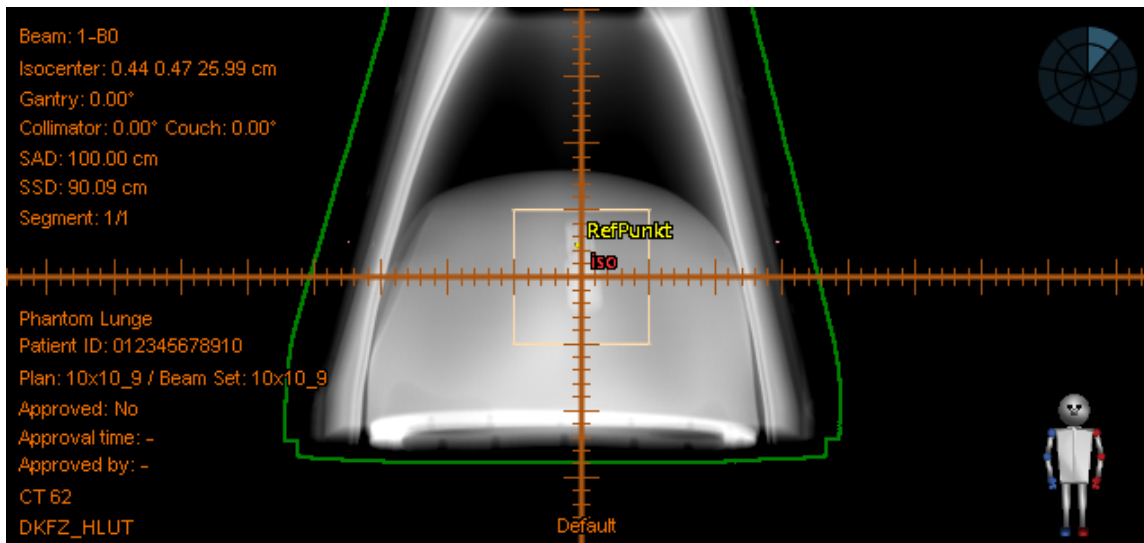

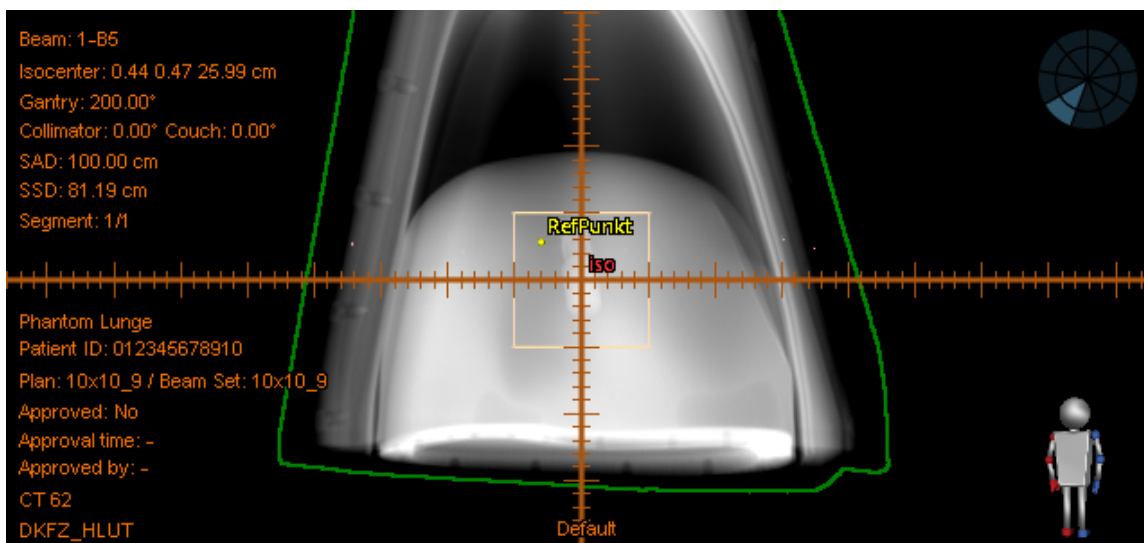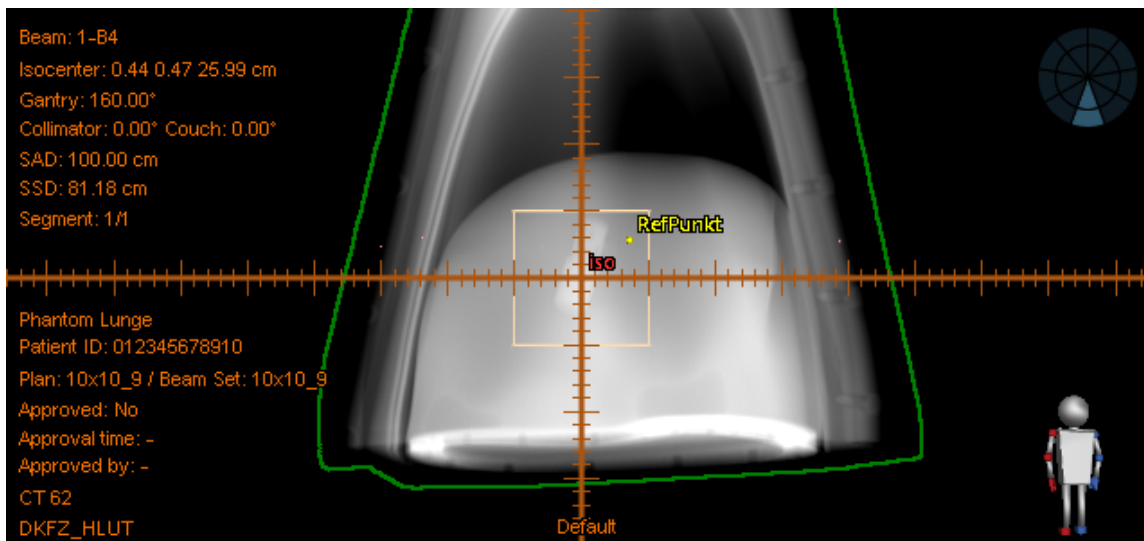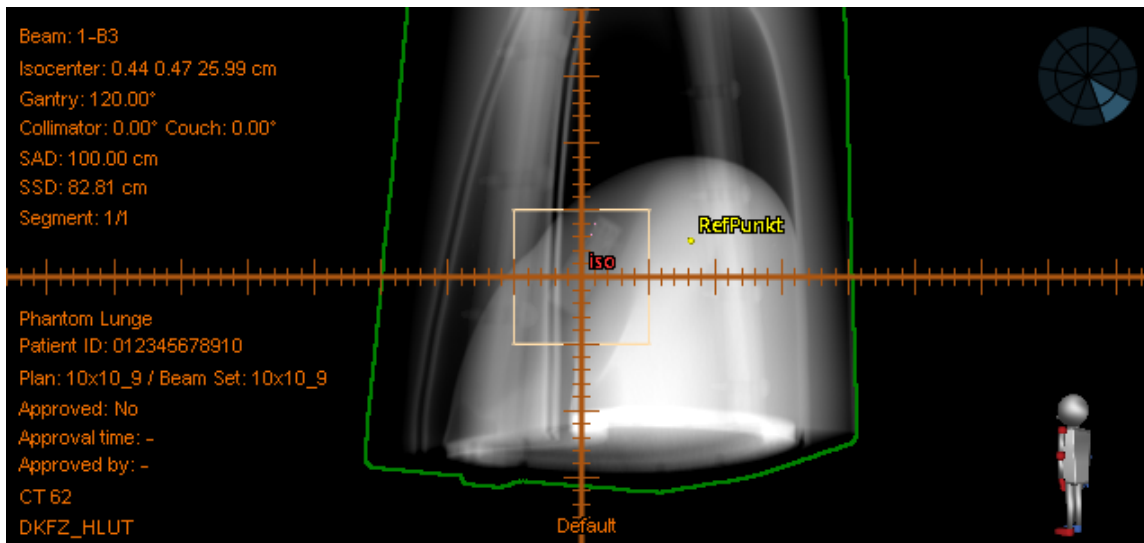

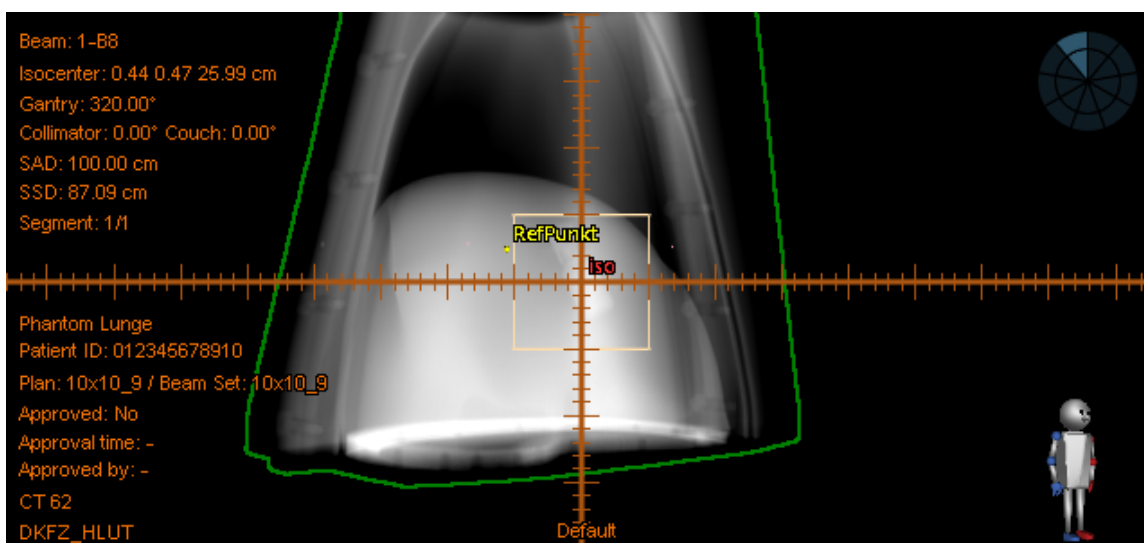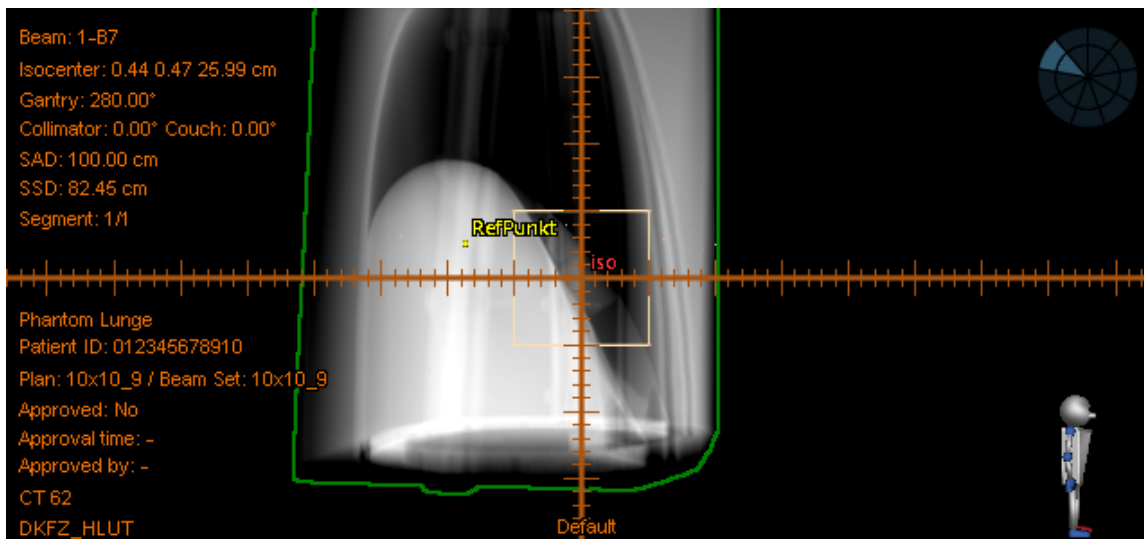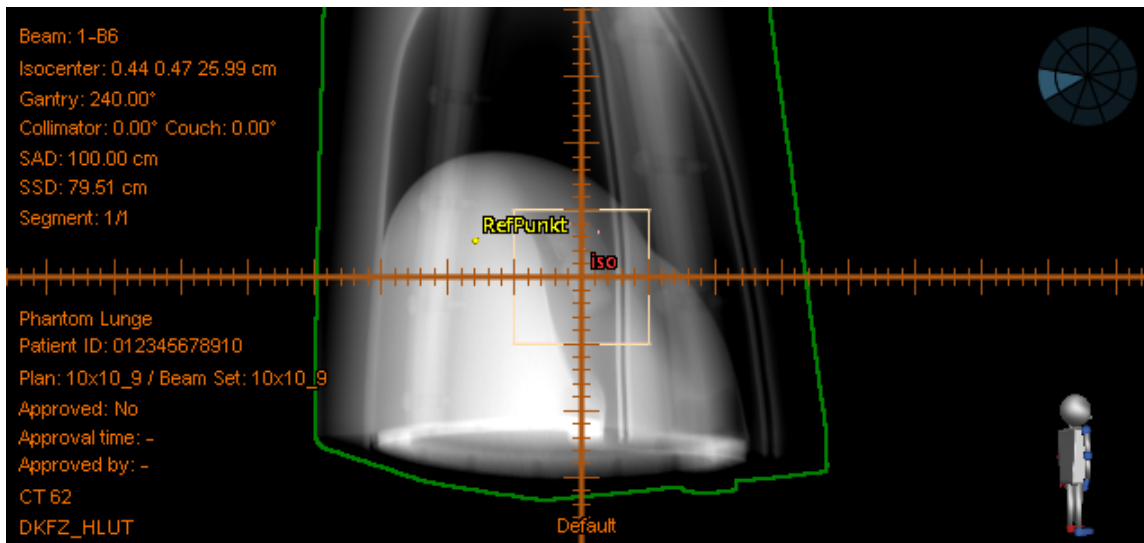



Patient name Phantom Lunge
Patient ID 012345678910
Treatment plan name 10x10_9

Plan last save time 13 Aug 2015, 09:21:49 (hr:min:sec)
Report creation time 13 Aug 2015, 16:35:24 (hr:min:sec)
Plan and structure set approved No
Plan approved by -
Plan approval time -

Import log

13 May 2015, 12:59:57 (hr:min:sec)	AD\meduser	Starting import. RayStation version 4.5.0.19
13 May 2015, 13:00:03 (hr:min:sec)	AD\meduser	DICOM import succeeded
13 May 2015, 13:11:13 (hr:min:sec)	AD\meduser	Starting import. RayStation version 4.5.0.19
13 May 2015, 13:12:10 (hr:min:sec)	AD\meduser	DICOM import succeeded
13 May 2015, 13:12:10 (hr:min:sec)	AD\meduser	Please note the following warnings / assumptions:
13 May 2015, 13:12:10 (hr:min:sec)	AD\meduser	#1 Patient's Name differed from the current patient.

Name of current patient:
 Lunge^Phantom

Mismatching names from imported data:
 LUNGE^PHANTOM

10 Jun 2015, 15:17:56 (hr:min:sec)	AD\meduser	Starting import. RayStation version 4.5.0.19
10 Jun 2015, 15:18:01 (hr:min:sec)	AD\meduser	DICOM import succeeded
10 Jun 2015, 15:18:01 (hr:min:sec)	AD\meduser	Please note the following warnings / assumptions:
10 Jun 2015, 15:18:01 (hr:min:sec)	AD\meduser	#1 Patient's Name differed from the current patient.

Name of current patient:
 Lunge^Phantom

Mismatching names from imported data:
 LUNGE^PHANTOM

17 Jun 2015, 14:31:35 (hr:min:sec)	AD\meduser	Starting import. RayStation version 4.5.0.19
17 Jun 2015, 14:31:38 (hr:min:sec)	AD\meduser	DICOM import succeeded
17 Jun 2015, 15:21:43 (hr:min:sec)	AD\meduser	Starting import. RayStation version 4.5.0.19
17 Jun 2015, 15:22:28 (hr:min:sec)	AD\meduser	DICOM import succeeded
17 Jun 2015, 15:22:28 (hr:min:sec)	AD\meduser	Please note the following warnings / assumptions:
17 Jun 2015, 15:22:28 (hr:min:sec)	AD\meduser	#1 Patient's Name differed from the current patient.

Name of current patient:
 Lunge^Phantom

Mismatching names from imported data:
 LUNGE^PHANTOM

29 Jul 2015, 13:17:38 (hr:min:sec)	AD\meduser	Starting import. RayStation version 4.5.0.19
29 Jul 2015, 13:17:44 (hr:min:sec)	AD\meduser	DICOM import succeeded
29 Jul 2015, 15:42:36 (hr:min:sec)	AD\meduser	Starting import. RayStation version 4.5.0.19
29 Jul 2015, 15:42:49 (hr:min:sec)	AD\meduser	DICOM import succeeded
29 Jul 2015, 16:09:27 (hr:min:sec)	AD\meduser	Starting import. RayStation version 4.5.0.19
29 Jul 2015, 16:11:21 (hr:min:sec)	AD\meduser	DICOM import succeeded
29 Jul 2015, 16:11:21 (hr:min:sec)	AD\meduser	Please note the following warnings / assumptions:
29 Jul 2015, 16:11:21 (hr:min:sec)	AD\meduser	#1 Patient's Name differed from the current patient.

Name of current patient:
 Lunge^Phantom

Mismatching names from imported data:
 LUNGE^PHANTOM

05 Aug 2015, 08:28:48 (hr:min:sec)	AD\meduser	Starting import. RayStation version 4.5.0.19
05 Aug 2015, 08:28:50 (hr:min:sec)	AD\meduser	DICOM import succeeded
05 Aug 2015, 09:18:38 (hr:min:sec)	AD\meduser	Starting import. RayStation version 4.5.0.19
05 Aug 2015, 09:20:25 (hr:min:sec)	AD\meduser	DICOM import succeeded
05 Aug 2015, 09:20:25 (hr:min:sec)	AD\meduser	Please note the following warnings / assumptions:
05 Aug 2015, 09:20:25 (hr:min:sec)	AD\meduser	#1 Patient's Name differed from the current patient.

Name of current patient:
 Lunge^Phantom

Mismatching names from imported data:
 LUNGE^PHANTOM

05 Aug 2015, 14:24:49 (hr:min:sec)	AD\meduser	Starting import. RayStation version 4.5.0.19
05 Aug 2015, 14:25:01 (hr:min:sec)	AD\meduser	DICOM import succeeded
05 Aug 2015, 14:52:21 (hr:min:sec)	AD\meduser	Starting import. RayStation version 4.5.0.19
05 Aug 2015, 14:54:15 (hr:min:sec)	AD\meduser	DICOM import succeeded
05 Aug 2015, 14:54:15 (hr:min:sec)	AD\meduser	Please note the following warnings / assumptions:
05 Aug 2015, 14:54:15 (hr:min:sec)	AD\meduser	#1 Patient's Name differed from the current patient.

Name of current patient:
 Lunge^Phantom

Mismatching names from imported data:
 LUNGE^PHANTOM



Patient name Phantom Lunge
Patient ID 012345678910
Treatment plan name 10x10_9

Plan last save time 13 Aug 2015, 09:21:49 (hr:min:sec)
Report creation time 13 Aug 2015, 16:35:24 (hr:min:sec)
Plan and structure set approved No
Plan approved by -
Plan approval time -

13 Aug 2015, 09:01:35 (hr:min:sec) AD\meduser
13 Aug 2015, 09:01:39 (hr:min:sec) AD\meduser
13 Aug 2015, 09:01:39 (hr:min:sec) AD\meduser
13 Aug 2015, 09:01:39 (hr:min:sec) AD\meduser

Starting import. RayStation version 4.5.0.19
DICOM import succeeded
Please note the following warnings / assumptions:
#1 Patient's Name differed from the current patient.

Name of current patient:
Lunge^Phantom

Mismatching names from imported data:
LUNGE^PHANTOM